


VD-06 Digital Velocity Decoder

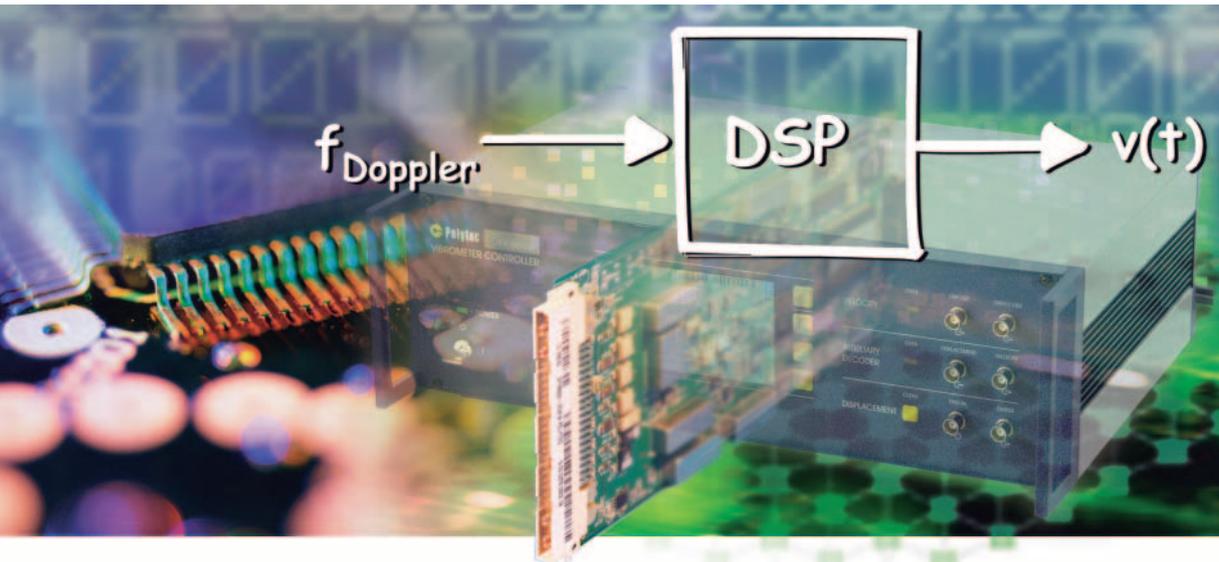

Modular Vibrometer System

- OFV-5000
Vibrometer Controller
– Velocity Decoders
– Displacement Decoders
- OFV-505/503
Standard Sensor Heads
- OFV-534
Compact Sensor Head
- OFV-551/552
Fiber Interferometers

Measuring vibration velocity

Signal processing is one of the most sensitive parts of any Laser Doppler Vibrometer system. The quality of dedicated velocity decoders available for the OFV-5000 Controller defines the accuracy, linearity, sensitivity and signal-to-noise ratio of the system.

Velocity Decoding in Laser Vibrometer Systems

Polytec Laser Doppler Vibrometers operate on the Doppler principle, measuring back-scattered laser light from a vibrating structure, to determine its vibrational velocity and displacement.

A vibrometer system comprises controller electronics and a non-contact standard-optic or fiber-optic sensor head. The controller provides signals and power for the sensor head, and processes the vibration signals. These are electronically converted by specially developed decoders within the controller to obtain velocity and displacement information about the test structure.

A vibrometer system based on the OFV-5000 Controller can measure vibration frequencies from 0 Hz up to 24 MHz. Velocities can be as small as 0.1 $\mu\text{m/s}$ up to 10 m/s. Different measurement ranges demand appropriate decoders. To meet this demand Polytec offers a range of analog and digital decoders with different characteristics. Two velocity decoders can be operated simultaneously.

The VD-06 Digital High-Precision Velocity Decoder

The VD-06 decoder uses digital signal processing (DSP) technology to enhance its performance. Vibration measurements can be taken with exceptional precision and resolution. The VD-06 can be installed side-by-side with the popular VD-02 decoder to give an optimum decoder combination.

Key features of the VD-06

- DSP-based, high precision velocity decoder
- Large bandwidth 0 to 350 kHz
- Outstanding linearity and frequency response
- 4 velocity ranges:
10 mm/s, 20 mm/s, 100 mm/s and
500 mm/s full scale peak
- Best resolution: 0.01 $\mu\text{m s}^{-1} / \sqrt{\text{Hz}}$


VD-06 Technical Data

Measurement Range	Full Scale Output (Peak)*	Resolution**	Signal Frequency Range	Max. Acceleration	Max. Linearity Error
mm/s/V	m/s	$\mu\text{m s}^{-1} / \sqrt{\text{Hz}}$	Hz	g	% of reading
1	0.01	0.01	0 – 20 k	128	0.1
2	0.02	0.05	0 – 350 k	4,500	0.1
10	0.10	0.05	0 – 350 k	22,000	0.1
50	0.50	0.06	0 – 350 k	110,000	0.1

* The full scale values correspond to the maximum output voltage of $10 V_{\text{peak}}$.

** The resolution is defined as the signal amplitude (rms) that produces 0 dB signal/noise ratio with 1 Hz spectral resolution at 50 % f_{max} .

VD-06 Range Diagram

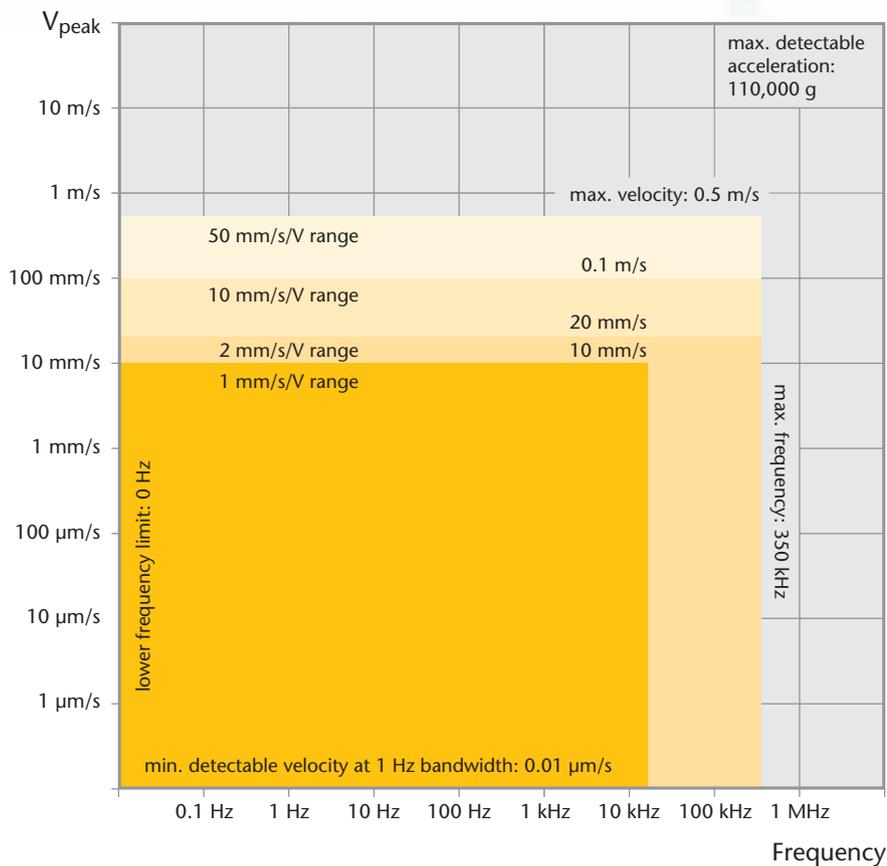

Data Acquisition

For PC-based data acquisition and processing, VibSoft software is recommended. VibSoft-1000 and VibSoft-1004 offer comprehensive software/hardware packages for dual and quad channel data acquisition from 0 to 1 MHz. VibSoft-80 and

VibSoft-84 are available for lower frequencies and cover dual and quad channels up to 80 kHz. For more information on selection and combination of signal decoders please see OFV-5000 and Decoder Guide data sheets, or contact your local vibrometer sales/application engineer.

Polytec GmbH (Germany)
 Polytec-Platz 1-7
 76337 Waldbronn
 Tel. +49 7243 604-0
 Fax +49 7243 69944
 info@polytec.de

Polytec France S.A.S.
 Bâtiment Orion – 1^{er} étage
 39, rue Louveau
 92320 Châtillon
 Tel. +33 1 496569-00
 Fax +33 1 57214068
 info@polytec.fr

Polytec Ltd. (Great Britain)
 Lambda House, Batford Mill
 Harpenden, Herts AL5 5BZ
 Tel. +44 1582 711670
 Fax +44 1582 712084
 info@polytec-ltd.co.uk

Polytec Japan
 Arena Tower, 13th floor
 3-1-9, Shinyokohama,
 Kohoku-ku, Yokohama-shi,
 Kanagawa, 222-0033
 Tel. +81 45 478-6980
 Fax +81 45 478-6981
 info@polytec.co.jp

Polytec, Inc. (USA)
 North American Headquarters
 16400 Bake Parkway
 Suites 150 & 200
 Irvine, CA 92618
 Tel. +1 949 943-3033
 Fax +1 949 679-0463
 info@polytec.com

Central Office
 1046 Baker Road
 Dexter, MI 48130
 Tel. +1 734 253-9428
 Fax +1 734 424-9304

East Coast Office
 25 South Street, Suite A
 Hopkinton, MA 01748
 Tel. +1 508 417-1040
 Fax +1 508 544-1225